

# Affinity™ 4 Birthing bed

Delivering Safety



Enhancing outcomes for  
patients and their caregivers:

**Hill-Rom**

# Our birthing bed is the perfect combination of safety and ease of use



## The only birthing bed you'll need

In labour and delivery, every second counts. You must be able to respond quickly and easily, especially in an emergency.

The Affinity 4 birthing bed delivers on the promise of a safe birthing bed that functions in the busiest environment with ease. No more hassle with a removable foot section. With the Affinity 4 **Stow and Go™** feature, the integrated foot section simply slides under the bed for birth and back into place afterward. Or, if you choose our ergonomically designed lift-off version, simply place the foot section on the built-in stand.



Ergonomic push handles and integrated permanent IV poles.



### Safety for patients and caregivers

- Foot section easily slides back into place
- **EasyGlide™** calf supports
- Four dual-locking casters to keep the bed securely in place
- Hill-Rom® high quality foam mattress for comfort and pressure redistribution
- Automatic pelvic tilt to reduce patient migration
- Clean perimeter around the bed with no protruding pieces
- Optional ergonomically designed lift-off section

### Ease of use

- No-lift, **Stow and Go™** foot section
- Battery backup
- Zero transfer gap
- Easy-to-use gas-assist Trend mechanism
- **OneStep™** siderail release

### New Options

- Ergonomic push handles with integrated IV pole cradle
- Integrated line managers
- Integrated permanent IV poles

# The Affinity 4 bed is designed to make every movement fast and easy

## Less time, less hassle

In labour and delivery, staff time is at a premium. No one can spare even a few minutes to search for equipment, look for an outlet or use complicated controls. The Affinity 4 bed is designed to reduce the hassles you face every day.



- The **Stow and Go** foot section slides and stores under the bed – no lifting required. The foot section easily slides back into place without false latching or interference with linens.
- With **battery backup**, you can use all the bed's features, including changing height and position, when you're away from an outlet, reducing the time from "decision to incision" in critical situations.
- The gas-assist Trend mechanism allows head elevation to immediately lower to 8 degrees.
- **OneStep** siderail release provides zero transfer gap for safer and easier patient transfers.
- The bed's streamlined design, with all the equipment stored completely within the perimeter (instead of protruding out), protects you as you move quickly around your patient.

# The Affinity 4 bed brings dozens of safety features to the bedside

## Safety and convenience

We have designed every detail of the Affinity 4 birthing bed to help make the birthing experience safer for you and your patient.

- The foot section locks firmly and easily in position.
- **EasyGlide** calf supports swing into the correct position and lock with a single lever, staying firmly in place.
- Trend mechanism, with improved springs and seals, makes lowering and lifting an easy, one-person job.
- Inflatable air pillows in back and seat decrease the need to reposition the patient for delivery.
- Four dual-lock casters keep the bed in place.
- Hill-Rom Comfortline™ foam mattress options for pressure redistribution and superior comfort, especially on the vulnerable seat area.
- Automatic pelvic tilt and optional v-cut mattress provide better access to patient for birth.



An optional lift-off foot section is available. Lightweight and easy to handle, the foot section can be stored on a built-in stand.



**EasyGlide** calf supports are larger, with a strong release lever that can be actuated from any angle with one hand.



Patient-friendly grips and siderails make it easier for patients to get up and moving.



Automatic pelvic tilt to 15 degrees helps keep the patient in place during labour and improves positioning during birth.

# Progressive Labour


with the Affinity 4 Birthing Bed



*Throne*



*Fours*

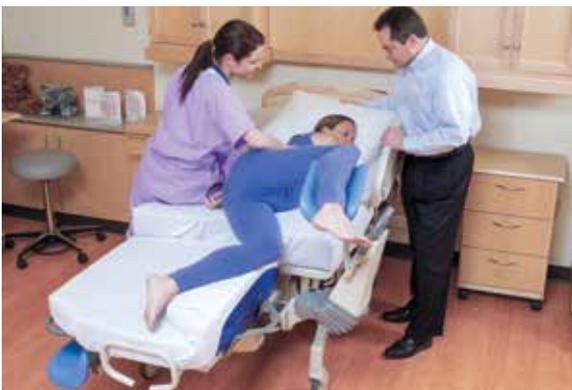

*Side*



*Squat*



*Tug-of-war*



*Fowlers*

## Technical Specifications

<b>Product Identification</b>	P3700 Affinity 4 Birthing Bed	<b>Critical Angles</b>	
<b>Bed Length</b>		Maximum head elevation	63°
Full	229 cm (90")	Maximum seat elevation	15°
<b>Bed Width</b>		Maximum Trend-like positioning	8°
Siderails up	99 cm (39")	<b>Maximum Lift Capacity</b>	
Siderails down	91 cm (36")	Bed	227 kg (500 lbs)
<b>Bed Height</b>		Foot section	181 kg (400 lbs)
Low (with mattress)	58 cm (23")	Caster size	15 cm (6")/20 cm (8")
High (with mattress)	99 cm (39")	Bed weight	218 kg (481 lbs)
Maximum height of seat section in Trend-like position	102 cm (40")		
<b>Mattress</b>			
Length	198 cm (78")		
Width	89 cm (35")		
Thickness (head section)	12.7 cm (5")		
Thickness (foot section)	8.9 cm (3.5")		

### Standard Features

- Battery Backup
- Standard 15 cm (6") Caster
- Automatic Night Light
- Lift-off or Stow and Go Foot Section
- Automatic Seat Tilt
- V-cut or Straight Edge Surface
- 4 Tente Locking Casters
- Blow-Molded Headboard
- Fluid Basin
- Central Brake and Steer
- Instant CPR Release
- Lockout Control and Indicator
- Instant Labour Grips
- IV Pole
- OneStep Siderail Release
- Trend-like Positioning
- Foot Section Removal
- Foot Supports

### Accessories

- Anaesthesia Screen
- ISS System
- Arm Board
- Labour Bar
- Calf Supports – Attached
- O<sub>2</sub> Tank Holder
- Calf Supports – Telescoping
- Pendant Control
- Disposable Drapes
- Permanent IV Pole
- Foley Hook Kit
- Sheets
- Full Leg Supports – Attached
- Slip Cover – Foot Section
- Full Leg Supports – Telescoping
- Slip Cover – Head Section

### Options

- Air Support System
- Variety of Wood Headboards
- Voltage – 120A or 230A
- Sidecom™ Communication Controls\*
- Comfortline Foam Mattress
- Ergonomic push handles with integrated IV pole cradle
- Integrated line managers
- Integrated permanent IV poles

Hill-Rom is a leading global medical technology company with more than 6,500 employees in over 100 countries. We partner with health care providers by focusing on patient care solutions that improve clinical and economic outcomes in five core areas: **Advancing Mobility, Wound Care and Prevention, Clinical Workflow, Surgical Safety and Efficiency, and Respiratory Health.** Hill-Rom people, programs, and product brands work towards one mission: **Every day, around the world, we enhance outcomes for our patients and their caregivers.**



Hill-Rom is committed to the concept that mobilising people early and supporting their independence, improves life - in the hospital and at home.

This document is destined solely for use by healthcare professionals. Medical devices shown in this brochure are intended for use with patients in departments of healthcare establishments.

**CE** These medical devices bear the CE mark and comply with all applicable regulations. Hill-Rom recommends you carefully read the detailed instructions for safe and proper use included in the documents accompanying the medical devices. The personnel of the healthcare establishments are responsible for the proper use and maintenance of these medical devices.

Medical beds manufacturer:

Hill-Rom SAS - B.P. 14 - Z.I du Talhouët - 56330 Pluvigner - France (also European Representative)

Völker GmbH - Wullener Feld 79 - 58454 Witten - Germany

Hill-Rom Holdings, Inc. - 1069 State Route 46 East - Batesville, IN 47006 - USA

Lifts manufacturer: Liko AB - Nedre vägen 100 - 975 92 Luleå - Sweden

Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.

©2014 Hill-Rom Services, Inc. ALL RIGHTS RESERVED.

Doc. No: 5EN137303-02, 03 September 2014

Not all products/options are available in all countries.

For further information about this product or a service, please contact your local Hill-Rom representative or visit our webpage:

[www.hill-rom.com](http://www.hill-rom.com)

Enhancing outcomes for  
patients and their caregivers:

**Hill-Rom**