


System Configuration

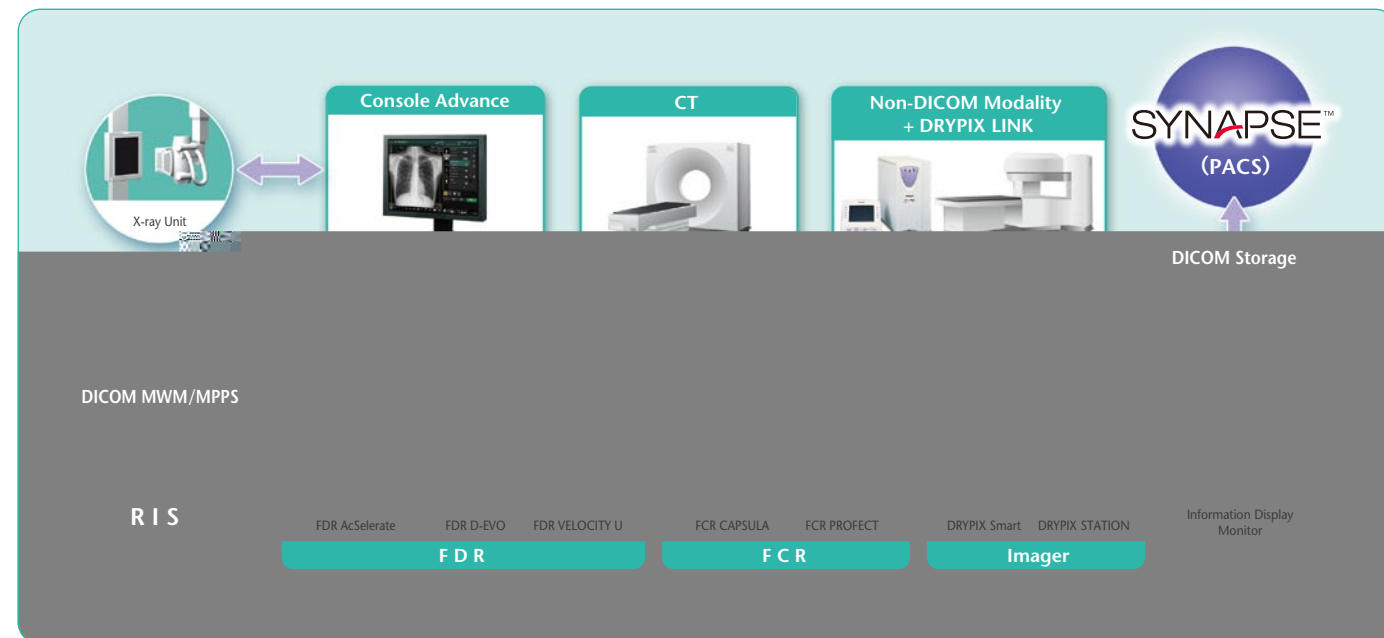


FUJI MEDICAL DRY LASER IMAGER

DRYPIX

Smart


Model: DRYPIX 6000

Highly efficient dry imager quickly offering excellent quality images for wider purposes

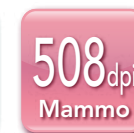

DRYPIX Smart Specifications

Standard Components	Fuji MEDICAL Dry Laser Imager DRYPIX Smart (Model: DRYPIX 6000)
Recording method	Laser exposure thermal development system
Applicable film	Fuji Medical Dry Imaging Film DI-HL (blue base) 35 × 43 cm (14" × 17") 35 × 35 cm (14" × 14") 26 × 36 cm (10" × 14") 25 × 30 cm (10" × 12") 20 × 25 cm (8" × 10") DI-ML (blue base) 26 × 36 cm (10" × 14") 25 × 30 cm (10" × 12") 20 × 25 cm (8" × 10")
Film loading	Daylight film loading
Film Tray	2 trays (5 sizes of film are available by changing film trays)
Processing capacity	Approx. 80 sheets/hour 35 × 43 cm (14" × 17")
Pixel size	50 μm (508 dpi)/100 μm (254 dpi)
Recording gradation	14 bits
Image memory	1GB
Density adjustment	Automatic
Input channels	DICOM network input × 1 channel only
Dimensions (W × D × H)	610 × 630 × 893 mm (24" × 25" × 35")
Weight	104 kg (229.3 lbs.)
Power Supply Conditions	Input voltage AC100 - 240V/ Single phase Frequency 50 - 60Hz
Environmental Conditions	Operating Conditions: • Temperature: 15 - 30°C • Humidity: 40 - 70%RH (at 15°C) to 15 - 70%RH (at 30°C) (No dew condensation)

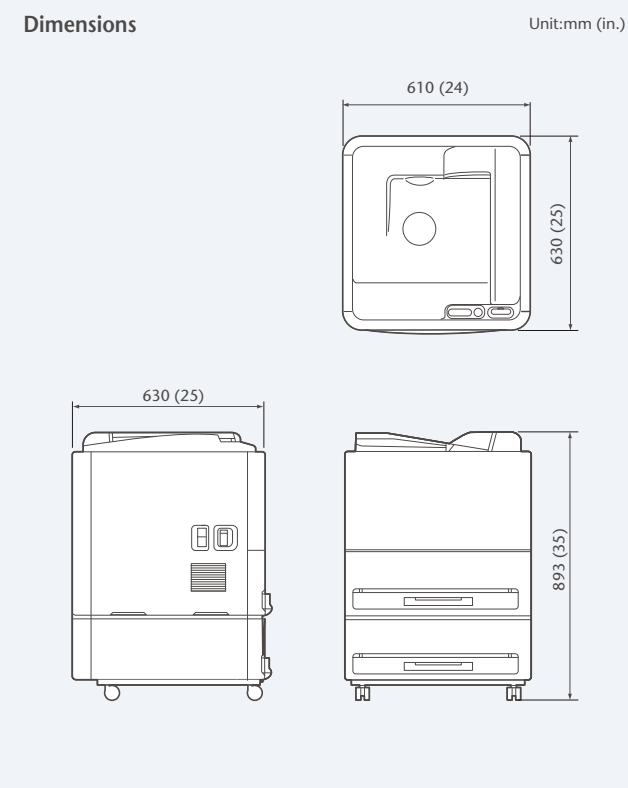

Specifications are subject to change without notice.
 All brand names or trademarks are the property of their respective owners.
 In some countries, regulatory approval may be required to import medical devices.
 For the availability of these products, please contact your local sales representatives.



CLASSIFIED
 UL US
 Medical Dry Imager
 WITH RESPECT TO ELECTRIC SHOCK, FIRE
 AND MECHANICAL HAZARDS ONLY
 IN ACCORDANCE WITH UL 60601-1/CAN/CSA C22.2 NO.6011
 EC 0050111 230K

CE
 European Authorized Representative:
 FUJIFILM Europe GmbH
 Hoosenstrasse 31, 40549 Duesseldorf, Germany

The most advanced DRYPIX has arrived, assisting smooth diagnoses

DRYPIX Smart, backed by Fujifilm's extensive experience in dry imaging, always delivers superior quality images to satisfy various needs of multi-department hospitals. Despite its compact size, enabling use anywhere in a medical facility, throughput is extremely high with no compromise on image quality.

Compact and highly efficient

High throughput

DRYPIX Smart boasts a world-class high throughput speed of 80 sheets per hour with 14" x 17" film. It will help reduce the patient's waiting time and greatly increase the efficiency of examination workflow.


Two trays to achieve more versatility

The DRYPIX Smart accommodates multiple film sizes. It is equipped with two universal film trays which enable printing on two different film sizes at the same time.



Fuji Medical Dry Imaging Film

The high quality DI-HL and DI-ML films contribute to producing clear images on the DRYPIX Smart. These films have a neutral color tone that produce images comparable to those made by wet processing.

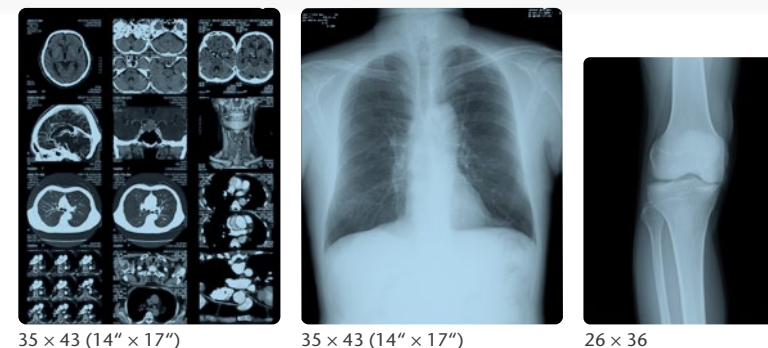

35 x 43 (14" x 17")

35 x 43 (14" x 17")

26 x 36

2 trays

Throughput
80 sheets/hr.

*with 14" x 17" film

Applicable for
mammography
(508 dpi)

High quality images for more versatility

High resolution and high maximum density


Offering high resolution of 508 dpi and a maximum density of 4.0*, the DRYPIX Smart is ideal for mammography which requires high definition images.

*When the DI-ML film is used.

Image processing engine which provides high-quality images


Advanced Variable Response (A-VR) Spline Interpolation Fujifilm's A-VR automatically detects and distinguishes between image data and alphanumeric characters, ensuring clear, sharp alphanumeric even when noisy images require smooth interpolation of image data. Benefits include easier, faster and more accurate diagnosis.

(A-VR) Flowchart

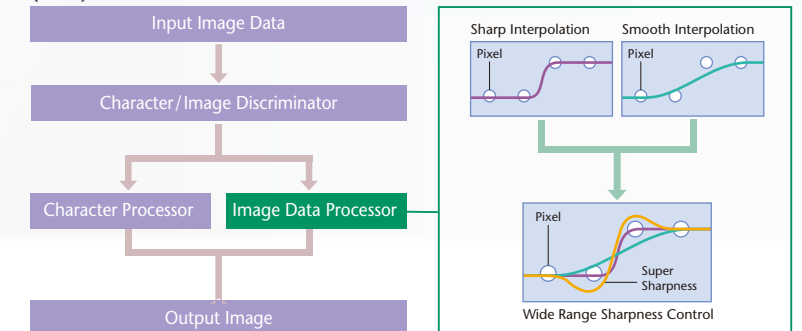

Quality Control

DRYPIX Smart prints a 24-step grayscale pattern to film, and then measures its density. This feedback system allows precise and subtle image adjustments (FDC: Auto Film Density Correction) to be made. Several kinds of test pattern images for the QC of mammograms are incorporated into DRYPIX Smart.

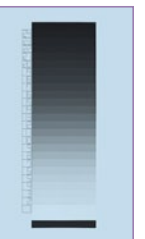

SAR (Smooth Curve Arranging)

Smooth Curve Arranging (SAR) on DRYPIX not only offers the most suitable image tones for modalities such as CT and MRI, but also allows adjustment of the tones to best match the diagnostic needs of individual patients. What's more, LUT also carries information on a wide range of modalities from different manufacturers to enable precise matching of image tone to specific modality.